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
# Naipo

## User Manual Shiatsu Massage Cushion with Heat for Lower Back

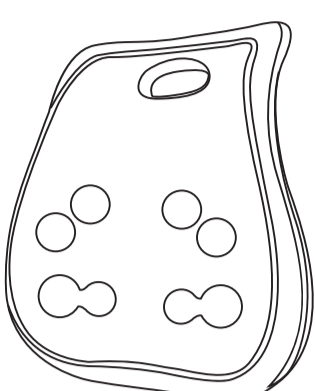Model: MGBK-136D

## Contents

English	01-03
Deutsch	04-06
Français	07-09
Español	10-12
Italiano	13-15
日本語	15-18

### Introduction



Thank you for purchasing the Naipo Shiatsu Massage Cushion with Heat for Lower Back. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

### Package Contents

- 1 x Shiatsu Massage Cushion
- 1 x User Manual
- 1 x Car Charger
- 1 x Charger

### Product Diagram

- Handle
- Massage Nodes
- Zipper
- Power Button
- DC Connector
- Power Plug
- Cable
- Appliance Plug
- Cigarette lighter
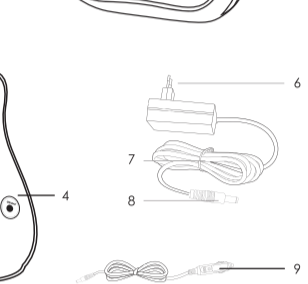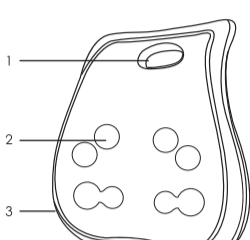
### Specifications

Dimensions (W*H*D)	36*14.5*33.5 cm
Weight	2 KG
Rated Output	30 W
Rated Voltage	DC 12V=2.5A

### Getting Started

- Plug the power adapter into the 220V / 50Hz outlet.
- Press the power button to activate the massage and heating function. Press it again to turn off the heating function. The selection of heating is cycles.
- Long press the power button for 3 seconds to turn off the device.
- The device will change the massage direction automatically per 1 min. (Outward / Inward)
- The device will stop automatically after 20 mins.
- Please turn off the device and unplug the power adapter from the outlet.

### Q&A

Problem	Cause	Remedy
Unable to start the device.	The cord is not connected properly.	Please push it harder until you heard the sound of "click".
The device doesn't work.	The usage time is too long.	Please stop using the product and unplug it from its power source for at least 1 hour.

### Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us at the address below that corresponds with your region. Please include your Amazon order number and product model number.

Website: www.naipocare.com  
Email: support@naipocare.com

\*Please note: Naipo can only provide after sale service for products purchased directly from Naipo. If you have purchased from a different seller, please contact them directly for service or warranty issues.

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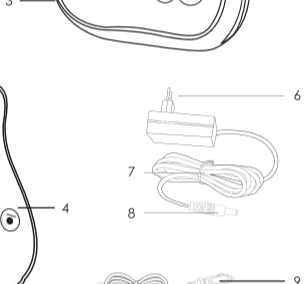
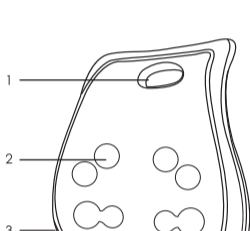
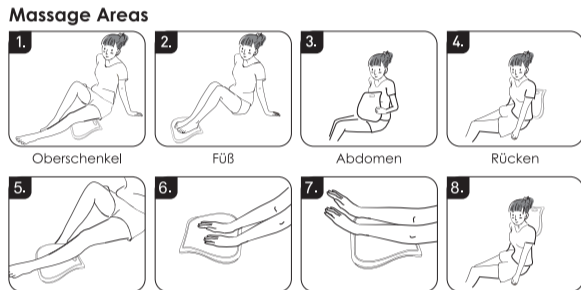
Website: www.czadeds.com  
Email-Adresse: support@czadeds.com

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