- | | |
|---|---|
| <p>(D) Maniküre/Pediküre-Set
Gebrauchsanweisung2</p> | <p>(RUS) Набор для маникюра и педикюра
Инструкция по применению ...59</p> |
| <p>(GB) Manicure/Pedicure set
Instruction for Use12</p> | <p>(PL) Zestaw do manicure/pedicure
Instrukcja obsługi.....70</p> |
| <p>(F) Ensemble de manucure/ pédicure
Mode d'emploi.....21</p> | <p>(NL) Manicure-/pedicureset
Gebruiksaanwijzing80</p> |
| <p>(E) Set de manicura/pedicura
Instrucciones de uso31</p> | <p>(P) Conjunto de manicure/ pedicure
Instruções de uso90</p> |
| <p>(I) Set per manicure/pedicure
Istruzioni per l'uso41</p> | <p>(GR) Σετ μανικιούρ/πεντικιούρ
Οδηγίες χρήσεως.....100</p> |
| <p>(TR) Manikür/Pedikür-Seti
Kullanma Talimatı.....50</p> | |

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Germany

ENGLISH

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Contents

1. Important notes	14	3.2 Attachments.....	17
2. Preparation	15	4. Cleaning/care	18
2.1 Information concerning the instrument	15	5. Disposal.....	19
2.2 Power supply	16	6. Replacement parts and wearing parts	19
3. Application.....	16	7. Warranty/service.....	20
3.1 General.....	16		

Dear customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, gentle therapy, blood pressure/diagnosis, weight, massage, beauty and air. Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

Items included in the package

- Manicure and pedicure device
- 7 high-quality attachments made from sapphire and felt
- Protective cap for nail dust
- Adapter
- Storage pouch
- These instructions for use

Warning

- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- Do not use the device if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- Repairs must only be carried out by Customer Services or authorised suppliers.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

1. Important notes

IMPORTANT:

- When in use, there may be an intense heat build-up on the skin or nail.
 - Check results regularly during treatment. This is particularly advisable in the case of diabetics, as they are less sensitive to pain and are liable to become injured more easily. Attachments 1, 2, 3, 4 and 5 that are provided are basically suitable for diabetics as they are only used for treatment of areas rather than specific points. Attachments 6 and 7 are, however, less suitable for diabetics. Always use with extreme care. If in doubt, ask your doctor.
 - This instrument must be used only for the purpose for which it has been developed, and in the manner described in the instructions for use. Any improper use can be dangerous.
 - In case of prolonged intensive use (e.g. milling of dead skin from your feet) the device might get very hot. In order to prevent skin burns in these situations, there should be longer breaks between each application. For your own safety, constantly monitor heat build up in the device. This particularly applies to people who are insensitive to heat.
 - The device is only for private use and is not intended for medical or commercial purposes.
 - For reasons of hygiene, the device should only be used on one person.
 - The manager is not responsible for damage or injury caused by improper or incorrect use.
 - Keep children away from the packaging materials. Risk of suffocation.
 - The appliance must be checked frequently for signs of wear or damage. If such signs are present, or if the appliance has been used improperly, it must be taken to either the manufacturer or the dealer before further use.
 - Switch off the instrument immediately in the event of defects or malfunctions.
 - Repairs may be carried out only by the customer service or authorized dealers.
 - On no account try to repair the instrument yourself!
 - Use the instrument only with the accessories provided.

- 
- Never operate the appliance without supervision, particularly when children are in the vicinity.
 - Do not use on small animals/animals.
 - Prevent all contact with water (except when cleaning with a slightly moistened cloth!). Water must never enter the inside of the appliance. Never immerse the appliance in water. On no account use the appliance in the bath, in a shower, in a swimming pool or over a water-filled wash basin. In the unlikely event that water gets into the case, disconnect the appliance immediately from the mains and contact your specialist electrical supplier or customer service.
 - Please ensure that no hair can get caught in the rotating attachments when using the device. Use a hair tie for your safety.
 - Keep the instrument away from sources of heat.
 - Do not use the instrument under blankets, pillows etc.
 - The instrument must not be operated continuously for longer than 20 minutes. After this time, wait an interval of at least 15 minutes to prevent the motor from overheating.
 - Ensure your hands are dry before handling either the instrument or the plug-in power unit.
 - Do not pull on the power cord or instrument to unplug the plug-in power unit from the socket outlet.
 - The plug-in power unit must be operated only from a power supply as specified on it.
 - The appliance should only be operated with the power unit provided.
 - After use, before cleaning and changing any accessory parts, the power unit should always be disconnected.

2. Preparation

2.1 Information concerning the instrument

This manicure and pedicure instrument is equipped with high-quality attachments. The basic version includes 7 attachments with long-life sapphire coating or polishing felt.

Together with the continuously adjustable speed regulation and clockwise/anti-clockwise rotation, professional foot and nail care at the standard of a pedicurist is possible. All the attachments can be easily integrated into the storage pouch of the MP 42.

The integrated LED light additionally creates ideal light conditions during application.

2.2 Power supply

To ensure adequate power at all times for good and powerful treatment of the feet and nails, this instrument is provided with a permanently attached power unit.

3. Application

3.1 General

- This device is intended for the treatment of hands (manicure) and feet (pedicure) only.
- Make sure that the device is switched off beforehand.
- Select the desired attachment part and carefully press it onto the axle of the device. To remove, pull apart the attachment and the device in opposing directions.
- Switch on the device by moving the sliding switch down (anti-clockwise rotation, LED lights up red) or up (clockwise rotation, LED lights up green).
- You can select the speed of the drive axle using the plus (+) and minus (-) buttons. Start each application at a low speed and then increase as required.
- Blue LEDs light up on the surface of the device depending on the selected speed.
- All filing and callus attachments are coated with a sapphire grain. This guarantees a long service life and ensures practically no wear.
- Do not press too hard and always carefully move the attachments towards the area to be treated.
- Move the device towards the area to be treated and apply slight pressure while slowly moving the device in circles.
- Please note that the callus attachments show reduced performance on soaked or wet skin. For this reason, we ask that you do not soak the foot in water before treatment.
- Do not remove all calluses to maintain natural skin protection.
- Ensure the axle can always rotate freely. The axle must not be permanently blocked, as otherwise the device will become too hot and may become damaged.
- Apply moisturiser to the treated areas after each use.

! Important:

Check the results during treatment. This is particularly important for diabetics, since their hands and feet are less sensitive.

3.2 Attachments

As these are attachments for use by professionals, we would ask you to handle them with care as they can cause injury when used incorrectly. Particular caution is advised at higher speeds.

The following attachments are provided with your appliance:

1 – Sapphire disk, fine

Filing and treatment of the nails, fine grained sapphire disk. The specific aspect of this sapphire disk is that only the inner disk rotates, while the outer support is stationary. This enables precise filing of the nails without the risk of burning the skin with the fast rotating disk.

2 – Sapphire disk, coarse

Filing and treatment of the nails, coarse grain sapphire disk. This attachment allows you to shorten even thick nails by filing. Since it quickly removes large areas of the nail, you should proceed with caution in this case also.

3 – Sapphire cone

Removal of dry skin or calluses from the soles of the feet and from the heels and for manicuring the nails.

4 – Felt cone

Smoothing and polishing the edge of the nail after filing as well as for cleaning the nail surface. Always polish in circular movements and do not allow the felt cone to rest in one place, otherwise friction could cause heat to develop.

5 – Sapphire callus rasp

For rapid removal of stubborn calluses and hard skin on the soles of feet and heels This attachment can be used on large areas.


6 – Cylindrical milling cutter

Grinding of lignified foot nail surfaces as well as smoothing of these. For this purpose, apply the milling cutter horizontally to the nail surface and remove the unwanted layer of nail using slow circular movements.

7 – Flame milling cutter

For removing ingrowing nails. Carefully guide the flame cutter up to the part to be treated and remove parts of the nail as necessary.

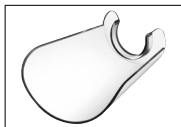8 – Protective cap for nail dust

This attachment ensures nail dust and skin cells are collected.

When filing, always work from the outside of the nail towards the tip!
Always check results regularly during treatment. As soon as the application becomes unpleasant, stop the treatment.

4. Cleaning/care

- Before cleaning the power unit, always unplug it from the socket outlet!
- Clean the instrument with a slightly moistened cloth. You can also moisten the cloth with mild soap suds to remove more stubborn dirt.
 - If necessary, for reasons of hygiene, the attachments can be cleaned with a cloth or brush moistened with alcohol. Dry them carefully afterwards.
 - Do not drop the instrument.

- Make sure that no water enters the instrument. But if this should happen, do not use the instrument again until it is completely dry.
- Do not use chemical detergent or scouring agents for cleaning.
- Never immerse the instrument or plug-in power unit in water or other liquids.

5. Disposal

In the interest of protecting the environment, the unit must not be thrown out with the household waste at the end of its service life.

Dispose of the unit at a suitable local collection or recycling point.

Observe the local regulations for material disposal.

Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).

If you have any queries, please contact the appropriate local authorities.

You can find out about collection points for your old units e.g. from the local municipal government, the local waste disposal company or your dealer.

6. Replacement parts and wearing parts

Replacement parts and wearing parts are available from the corresponding listed service address under the stated material number.

Designation	Item number and/or order number
Sapphire disk, fine	162.750
Sapphire disk, coarse	162.751
Sapphire cone	162.752
Felt cone	162.753
Sapphire callus rasp	162.754
Cylindrical milling cutter	162.755
Flame milling cutter	162.756
Protective cap for nail dust	163.755

7. Warranty/service

In case of a claim under the warranty please contact your local dealer or the local representation which is mentioned in the list “service international”.

In case of returning the unit please add a copy of your receipt and a short report of the defect.

The following warranty terms shall apply:

1. The warranty period for BEURER products is either 3 years or- if longer- the country specific warranty period from date of Purchase.
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Repairs (complete unit or parts of the unit) do not extend the warranty period.
3. The warranty shall not be valid for damages because of
 - a. improper treatment, e.g. nonobservance of the user instructions.
 - b. repairs or tampering by the customer or unauthorised third parties.
 - c. transport from the manufacturer to the consumer or during transport to the service centre.
 - d. The warranty shall not be valid for accessories which are subject to normal wear and tear (cuff, batteries etc.).
4. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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